


Sharpness. Consistency. Value.

SharpPoint Microsurgical Knives are designed to satisfy the needs of a full range of ophthalmic surgical procedures. They deliver outstanding precision and control resulting in accurate incisions, exceptional wound architecture and smooth, self-sealing wounds.

Unparalleled Sharpness

- Higher penetration capability than competitor offering
- Infinite edge guarantees industry lead technology that result in better control
- Creates unique wound architecture for optimal self-sealing incisions

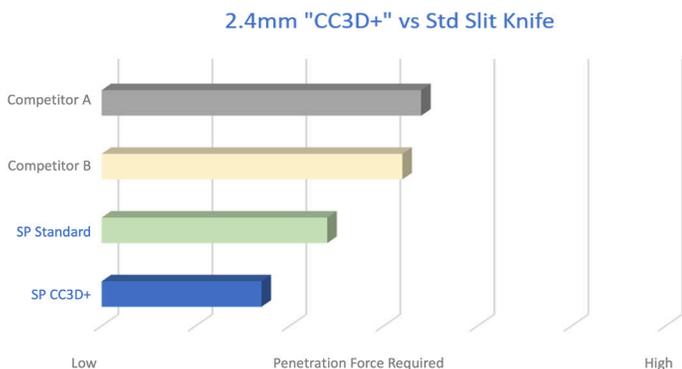

Proven Consistency


- Consistent, predictable and repeatable performance
- 50 years of guaranteed high quality

Comprehensive Value


- Comprehensive portfolio/partnership that includes Katena instruments SharpPoint suture, single-use cannulas and more!

* Data in chart demonstrates blade sharpness of the 4 leading manufacturers, utilizing our new in progress, Coefficient of Sharpness Test designed to verify blade performance

CC3D Plus Blade

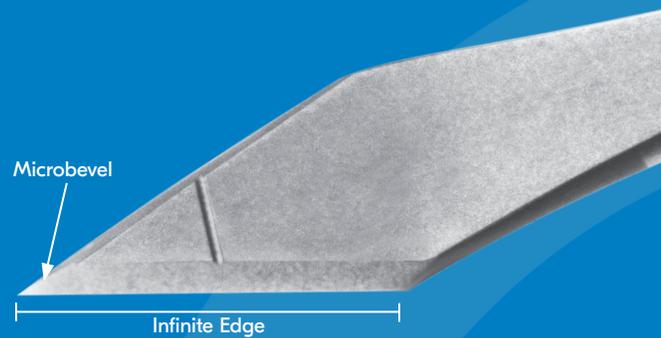

- Proprietary tip design and thinner blade prevents tenting during incision resulting in superior wound architecture
- Tip design adds vertical dimension, placing less stress on wound margins
- Microbevels provide stability and create a smooth, well-controlled, three-dimensional incision which places less stress on the wound margins, allows the wound to be more flexible during lens insertion with potentially few intra-corneal tears, and helps to promote rapid healing

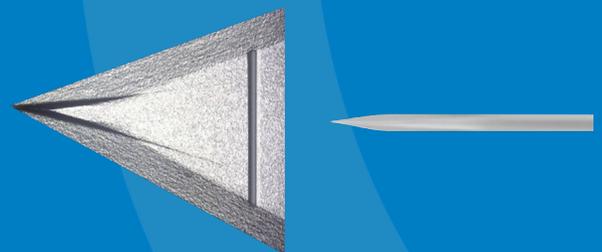

Top view, Magnification 100x

Side view, Magnification 100x

Infinite Edge Design = Sharper Blade!

- 33% less force needed to penetrate cornea tissue
- 22% thinner than our Standard slit knives equals even more control over incision
- Proprietary Infinite Edge™ Technology

Slit Knives



Phaco Incision Knives

Sharpoint angled precision slit knives penetrate easily while creating an accurate width incision. 6 per box.

Description	Single Bevel	Double Bevel
1.8mm Slit Angled		72-1831
2.2mm Slit Angled	72-2261	72-2231
2.4mm Slit Angled	72-2461	72-2431
2.5mm Slit Angled	72-2561	72-2531
2.65mm Slit Angled	72-2661	72-2631
2.75mm Slit Angled	72-2761	
2.8mm Slit Angled		72-2831
2.85mm Slit Angled	72-2861	
3.0mm Slit Angled	72-2961	72-2931
3.2mm Slit Angled	72-3261	72-3231
4.5mm Slit Angled	72-4561	

Specialty Knives



Clear Control 3D Plus

Sharpoint Clear Control 3D Plus slit knives are designed to facilitate instrument insertion through a multi-dimensional incision. Micro bevels are designed to add vertical dimension, placing less stress on wound margins. 6 per box.

Description	Double Bevel
1.6mm Angled	72-1650
1.8mm Angled	72-1850
2.0mm Angled	72-2050
2.2mm Angled	72-2250
2.4mm Angled	72-2450
2.5mm Angled	72-2550
2.65mm Angled	72-2650
2.75mm Angled	72-2750
2.85mm Angled	72-2850
3.0mm Angled	72-3050
3.2mm Angled	72-3250

Manufactured by Unique Technologies, Inc

Specialty Knives



UniBevel Clear Corneal Knives

Double-bevel slit knives create accurate width incisions easily with virtually no tissue drag or distortion. 6 per box.

Description	Double Bevel
2.2mm Angled	52-2242
2.4mm Angled	52-2442
2.5mm Angled	52-2542
2.65mm Angled	52-2642
2.75mm Angled	52-2742
2.85mm Angled	52-2842
3.0mm Angled	52-3042
3.2mm Angled	52-3242

Manufactured by Unique Technologies, Inc

Side Port



ClearPort Parenthesis Knives

ClearPort™ knives are designed to create consistent, accurate and self-sealing paracanthesis incisions. 6 per box.

Description	Angled	Straight
0.8mm, angled, full handle	78-2310	71-2301
1.0mm Angled	78-2020	71-2001
1.1mm Angled	78-2010	
1.4mm Angled	78-2019	

Side Port



Stab Knives

Sharpoint Stab Knives are heat-tempered to provide a strong and true point for easy penetration during initial stab incisions. 6 per box.

Description	Knife	Blade & Tip
15° Straight	72-1501	72-1502
22.5° Straight	72-2201	72-2202
30° Straight	72-3001	72-3002
45° Straight	72-4501	72-4502